



## The Assortment of Sterile Wound Dressings

### LEUKOMED®



Order Codes	Dressing Size	Wound Pad Size	Pieces / Box	Boxes / Case
7238014	2.8 in. x 2 in.	1.5 in. x 1.1 in.	50	10
7238001	3.1 in. x 4 in.	1.6 in. x 2.4 in.	50	8
7238002	3.1 in. x 6 in.	1.6 in. x 4.3 in.	50	8
7238003	4 in. x 7.8 in.	2 in. x 6.3 in.	50	8
7238004	4 in. x 9.8 in.	2 in. x 7.9 in.	50	8
7238005	4 in. x 11.8 in.	2 in. x 9.8 in.	50	8
7238006	4 in. x 13.8 in.	2 in. x 11.8 in.	50	8

### LEUKOMED® T




Order Codes	Dressing Size	Pieces / Box	Boxes / Case
7238106	2.8 in. x 2 in.	50	10
7238101	3.1 in. x 4 in.	50	10
7238102	4.3 in. x 5.5 in.	50	10

### LEUKOMED® T PLUS



Order Codes	Dressing Size	Wound Pad Size	Pieces / Box	Boxes / Case
7238212	2.8 in. x 2 in.	1.5 in. x 1.1 in.	50	10
7238201	3.1 in. x 4 in.	1.6 in. x 2.4 in.	50	8
7238202	3.1 in. x 6 in.	1.6 in. x 4.3 in.	50	8
7238203	3.9 in. x 9.8 in.	2 in. x 7.9 in.	50	6
7238204	3.9 in. x 11.8 in.	2 in. x 9.8 in.	50	5
7238205	3.9 in. x 13.8 in.	2 in. x 11.8 in.	50	4

## Perfect solutions for every wound

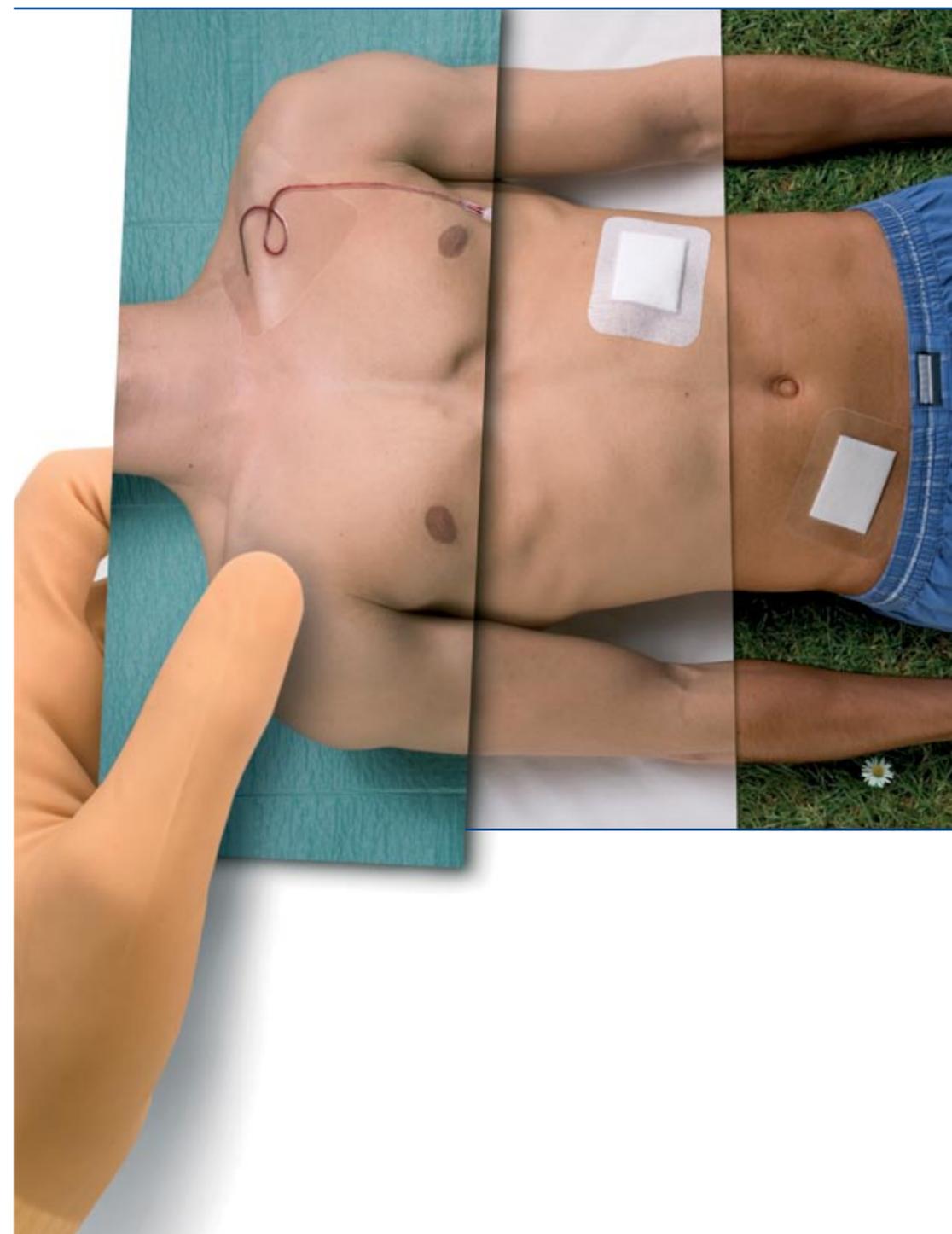

Sterile wound protection with  
**Leukomed®**  
**Leukomed® T**  
**Leukomed® T Plus**

53110 RN  
 L-LEUKOMED

# Leukomed®

## Sterile, Highly Absorbent Dressing for Heavily Secreting Wounds

Leukomed® is a dressing composed of an acrylic adhesive and non-woven substrate with a highly absorbent wound pad for secreting wounds.

### Indications

- Post surgical incision, especially on areas where swelling or movement is expected
- Moderate to heavy drainage wounds



Leukomed's wound pad has a **good absorption capacity, 12 times its weight**, allowing for long intervals of dressing changes and provides secure handling of wound fluid.



The **extra-soft** wound pad has a **non-adherent layer** that prevents it from sticking or interrupting the healing process.



The dressing is **permeable to air** thus reducing the risk of maceration.



# Leukomed® T

## Waterproof / Bacteria Barrier Film Dressing


Leukomed® T is a sterile, breathable, waterproof, transparent film dressing designed to cover catheter or wound sites while providing a barrier to external contaminants such as bacteria, blood and body fluids.

### Indications

- For use over catheter sites and tubes (peripheral and central)
- For use in clean, closed wound and skin protection

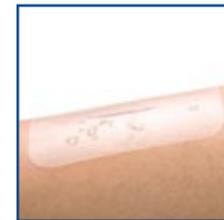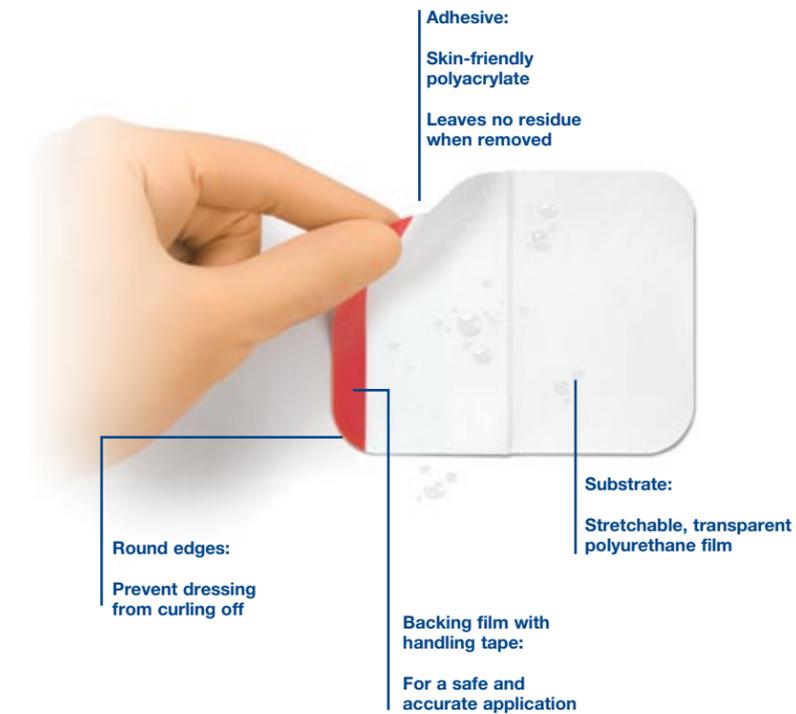

The **transparent polyurethane film** of Leukomed T is **waterproof and bacteria proof**, which protects the wound site even during showering or bathing.

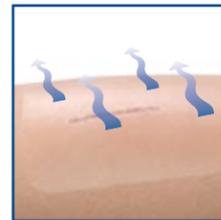

Leukomed T is **permeable to water vapor and oxygen**, reducing the risk of maceration.



**Easy-to use** patented delivery system.



# Leukomed® T Plus

## Highly Absorbent, Waterproof / Bacteria Barrier Film Dressing

Leukomed® T Plus is a sterile, breathable, waterproof, transparent bacteria barrier film dressing with a highly absorbent wound pad for secreting wounds.

### Indications

- Post-surgical dressing to protect clean, closed wounds
- Moderate to heavy wounds
- Waterproof to allow bathing

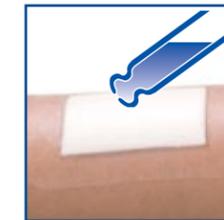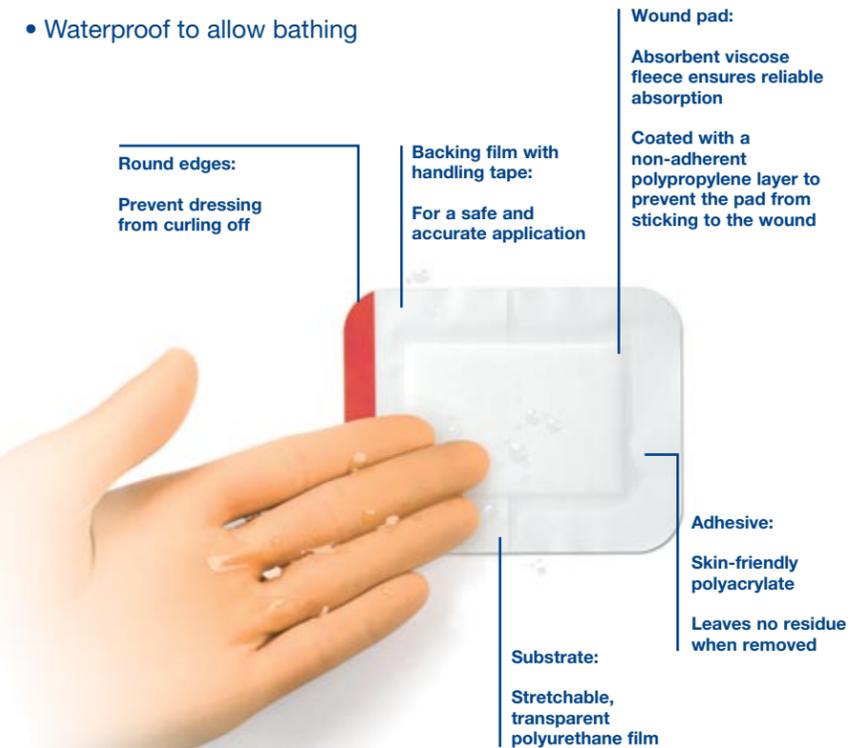

Leukomed T Plus' wound pad has a **good absorption capacity, 12 times its own weight**, allowing for long intervals of dressing changes and providing secure handling of wound fluid.



The **transparent polyurethane film** of Leukomed T Plus is **waterproof and bacteria proof**, which allows patients to shower and protect the site from outside contaminants.



The **extra-soft** wound pad has a **non-adherent layer** that prevents it from sticking or interrupting the healing wound.

### Application

1. Hold the dressing at the red handling tape. Then remove the upper backing paper on the non-red side.
2. Using the handling tape, the second backing paper can now be removed.
3. Fully press the transparent film onto the skin, especially making sure it is touching skin around the edges.
4. Take the red handling tape and remove the application film.

### Removal

1. Take one corner of film and pull – fully extending the stretch of the film. Painless and minimizes hair removal.